

INTRODUCING


MammoStar™

MammoStar™ biopsy site identifiers are the only organic long-term, ultrasound visible tissue markers available for the breast and axillary basin.

Combining ceramic, carbon and beta-glucan elements, MammoStar™ markers are an organic alternative to titanium or stainless steel markers.

The unique composition of MammoStar™ markers provides long-term visibility under Ultrasound*, X-ray and MRI.

*Remains visible for as long as 6 months.

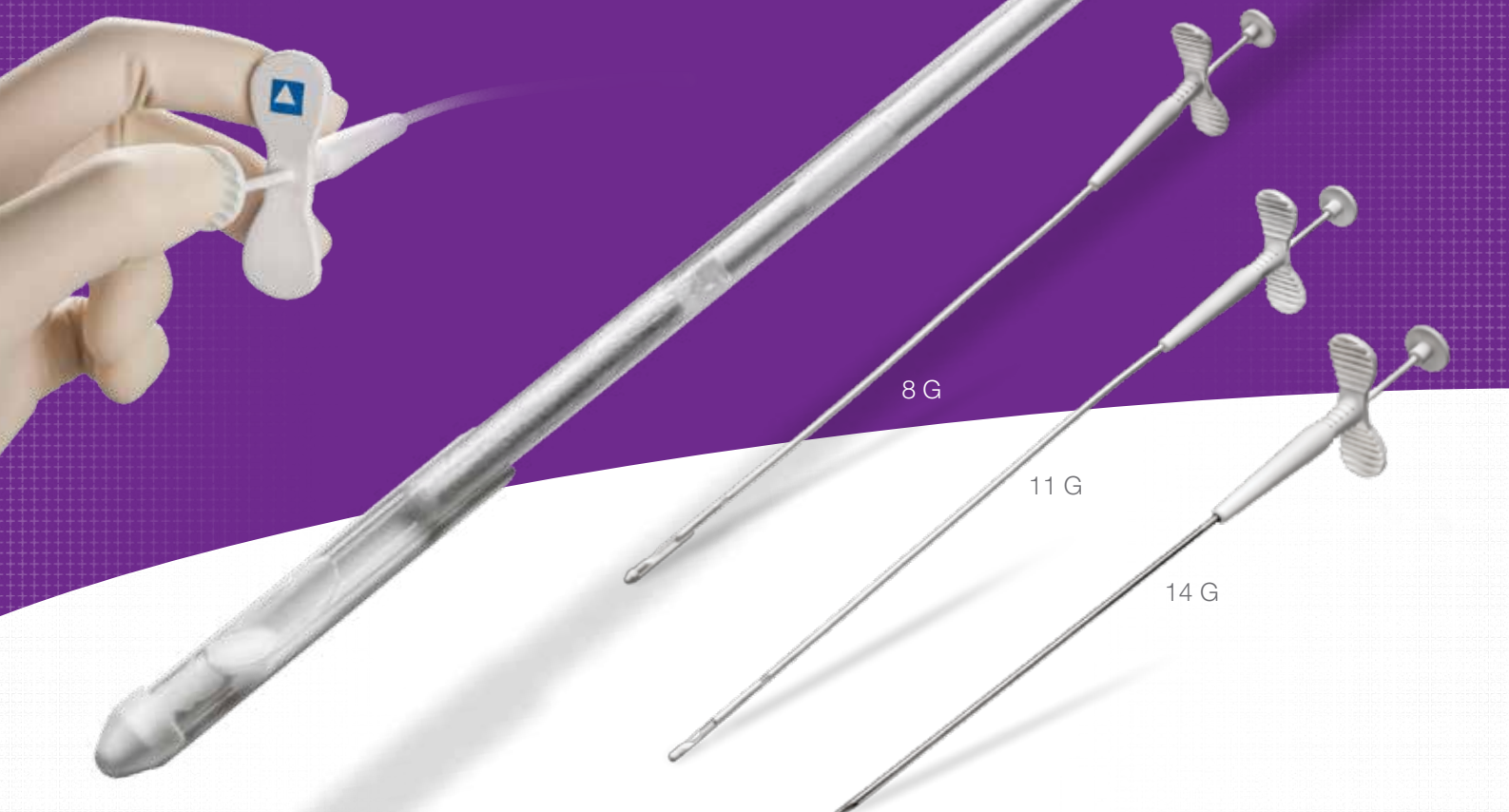

enduring. bright. distinct.

Mammotome® markers
biopsy site identifiers



Organic

MammoStar™ markers are comprised of three unique organic elements

Ceramic Marker

- + Chosen for biocompatibility, strength and unique imaging characteristics
- + Radiopaque
- + Non-ferromagnetic

Pyrolytic Carbon Coating

- + Enhanced ultrasound visibility
- + Carbon is found in the human body—helping to prevent rejection


Lyophilised Beta-Glucan Gel

- + Comprised of polysaccharide—a simple sugar
- + Absorbs surrounding fluid resulting in a gel-like material, suitable for ultrasound visibility
- + Pliable, Non-palpable
- + 5-year shelf life

Why Beta Glucan?

- + It is the carrier gel or suspending agent for the marker.
- + Beta glucans are sugars found in the cell walls of plants, such as oats.
- + All-organic. Not animal-derived.
- + Breaks down naturally over time.


Ultrasound Image

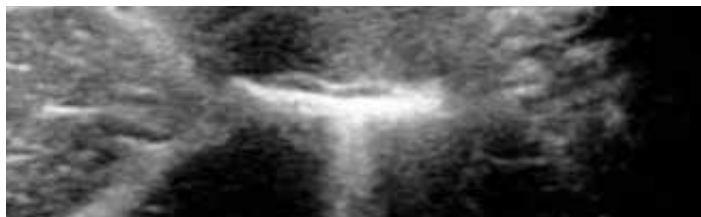

X-ray Image



MRI Image

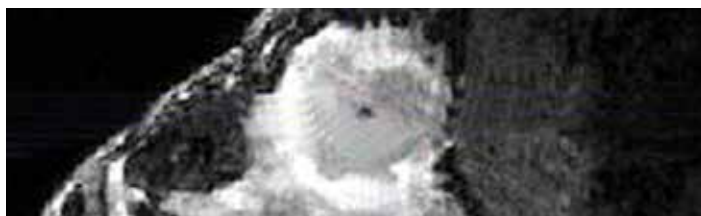

Order Placement

Simply choose your preferred option

+ Tel

Please contact Mammotome Customer Services at **0800 862 0587**.

+ Fax

Send your order via fax to Mammotome Customer Services at **0800 862 0585**.

+ Email

Send your order via email to Mammotome Customer Services at **orders.uk@mammotome.com**.

PRODUCT	CODE	DESCRIPTION	UNITS PER BOX
8 G	040525	MammoStar*	10
11 G	040526	MammoStar*	10
14 G	040527	MammoStar*	10

*Formerly known as BiomarC Plus from Carbon Medical Technologies

Mammotome
A Division of **DEVICOR** MEDICAL PRODUCTS

Devicor Medical UK Ltd. • 20-22 Bedford Row • London • WC1R 4JS • United Kingdom
Tel: 0800 8620587 • Fax: 0800 8620585 • Email: info.uk@mammotome.com • www.mammotome.com

Information for European countries outside of the UK by:
DEVICOR MEDICAL GERMANY GmbH • Südportal 3 • D-22848 Norderstedt
Tel: +49 (0) 40-59 35 59 10 • Fax: +49 (0) 40-500 98 940 • Email: info-europe@mammotome.com • www.mammotome.com

© 2012 Devicor Medical Products, Inc. All rights reserved. For complete product details, see Instructions for Use.
MDM 12-0022UK